


Distributed by:



# Guide to Reliable pH, Ion and Conductivity Measurements

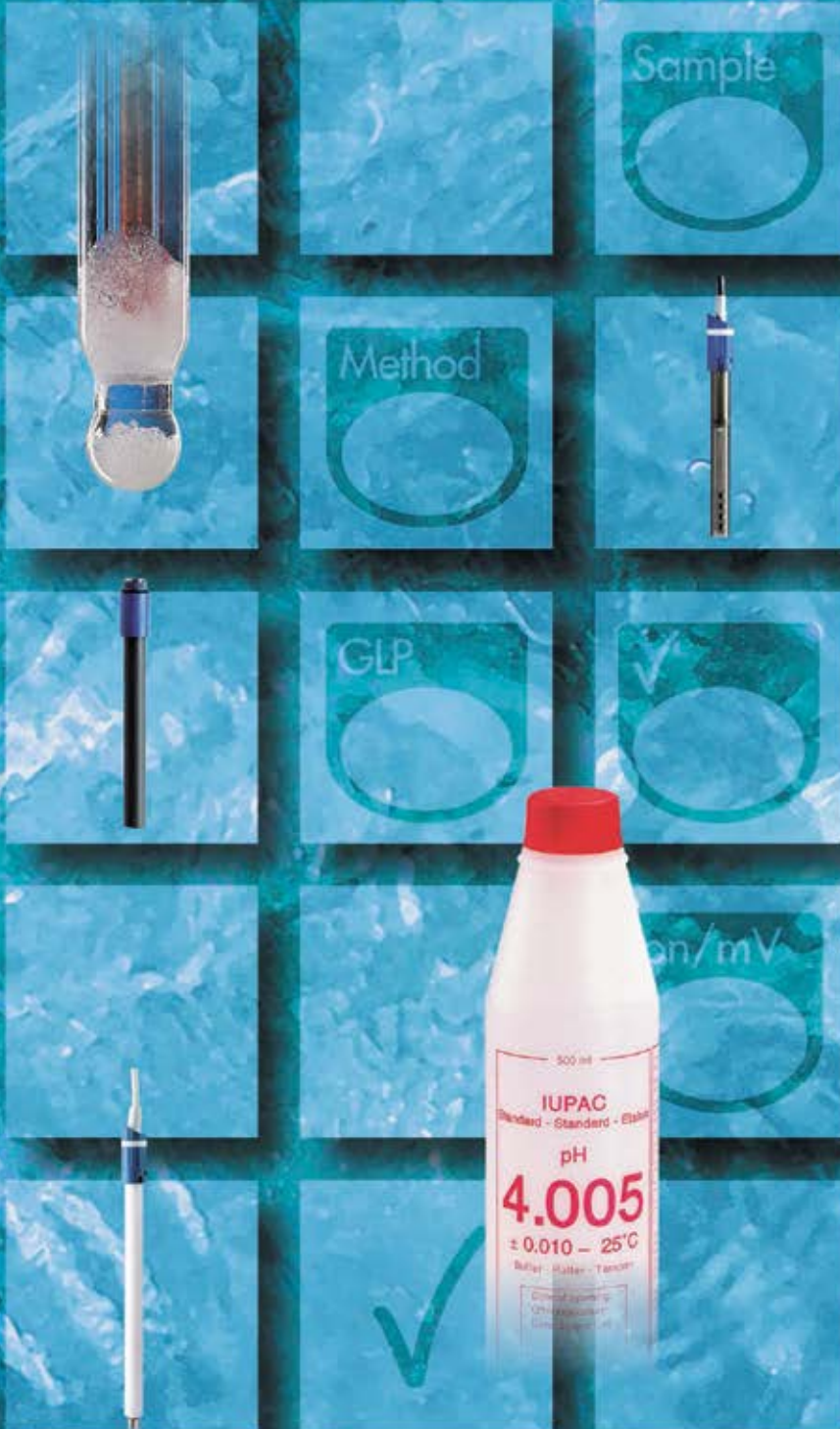

# The MeterLab® Concept

## Following Good Laboratory Practice is second nature with MeterLab

MeterLab from Radiometer Analytical incorporates all elements of the measuring chain - pH, ion and conductivity meters, electrodes and conductivity cells, solutions and sample handling - ensuring totally reliable measurements both in the laboratory and in the field. All elements are designed to make operation simple and error-free.



Radiometer Analytical offers a wide range of electrodes for every application and budget: combined pH, glass or reference electrodes featuring Red Rod or traditional technology, metal electrodes, ion-selective electrodes and conductivity cells. To choose the right pH electrode for your particular application, first refer to the table on pages 4 and 5.

For pH and conductivity measurements, Radiometer Analytical manufactures standards of the highest quality. Certain ranges are delivered with certificates of conformity and traceability to meet your accreditation and certification needs.

To keep your electrodes in excellent shape and extend their lifetime, Radiometer Analytical provides maintenance and filling solutions as well as necessary accessories.

Radiometer Analytical also offers a range of meters, sample stands and sample changers to help you ensure reliable pH, ion and conductivity measurements. For your full confidence, MeterLab laboratory meters now come with a free 5-year guarantee. To find out more, please ask for our separate MeterLab brochure or visit us at [www.hach-lange.com](http://www.hach-lange.com).



*The Right Electrode for your Application pages 4-5*



*Combined pH Electrodes pages 6-7*




*Reference Electrodes pages 8-9*



*Glass pH Electrodes page 10*



*Metal Electrodes page 11*

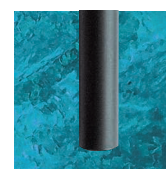

*Ion-selective Electrodes pages 12-13*



*Conductivity Cells pages 14-15*



*Solutions for Calibration and Maintenance pages 16-17-18*



*Accessories page 19*



## Red Rod Electrodes


**Radiometer Analytical's unique Red Rod Technology provides not only fast response time but also long-term accuracy and reproducibility.**

Red Rod electrodes feature:

**Reproducible measuring conditions** ensured by the special bayonet head. The electrode fits into place in one easy movement and is fixed safely and securely in the same position for every measurement.

**Wide temperature range** (0 to 100°C) thanks to fast temperature response of Red Rod Technology.


**Trouble-free liquid junction** as Red Rod encapsulation of reference elements means no silver ions and therefore less risk of clogging.

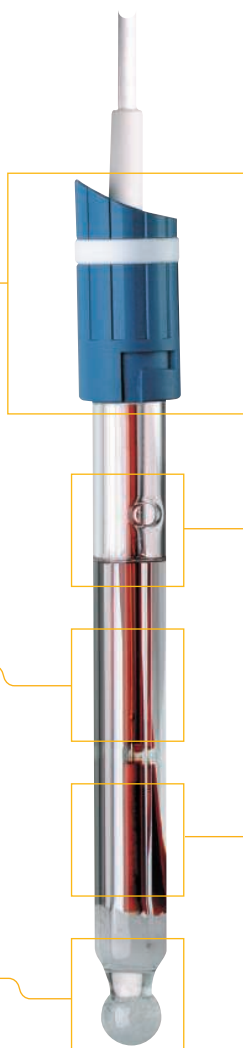

**Fast and accurate response** thanks to the symmetrical electrode chain formed by the internal and external Red Rod reference elements.

**Stable signal** as the Red Rod Technology protects the reference elements from light interference. Saturated KCl salt-bridge aids stable liquid junction potential.

### Your guarantee of reproducible results

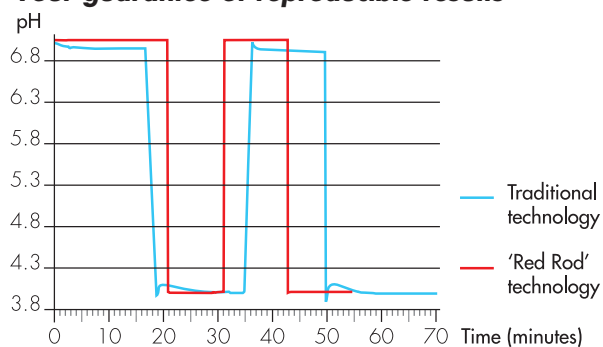

Red Rod combined pH electrode compared to traditional combined pH electrode

The symmetry of Red Rod electrodes means that the iso pH is the same as the zero pH which provides highly reproducible results even when temperature fluctuations occur.

The figure compares the measurements performed with a general-purpose Red Rod electrode and a traditional combined electrode. A MeterLab setup was used to measure certified standard pH solutions, pH 7.000 thermostated at 50°C and pH 4.005 thermostated at 25°C. The figure shows the fast response of the Red Rod electrode to changes in pH and temperature (curves offset to increase readability).

# The Right Electrode for your Application

Applications	Combined pH electrodes																																
	pHC2001	pHC2002	pHC2003	pHC2005	pHC2011	pHC2015	pHC2051	pHC2085	pHC2401	pHC2441	pHC2501	pHC2601	pHC2701	GK240TB	GK2401C	pHC3001	pHC3005	pHC3006	pHC3006L	pHC3011	pHC3031	pHC3081	pHC3085	pHC3105	pHC3131	pHC3185	pHC3359	pHC4000	pHC4001	pHC4006	XC161		
Alkaline solutions					✓	✓								✓							✓												
Aqueous solutions	✓	✓	✓	✓	✓	✓	✓	✓	✓			✓	✓	✓	✓	✓	✓	✓	✓	✓	✓		✓	✓	✓		✓	✓	✓	✓	✓	✓	✓
Blood												✓	✓																				
Education purposes																✓	✓																
Emulsions												✓																					
Fat/cream/cosmetics																						✓											
Field use				✓				✓									✓							✓	✓			✓					
Gel electrode																								✓	✓	✓	✓						
Hops/beer	✓	✓	✓	✓			✓	✓	✓				✓			✓	✓	✓					✓	✓									
Hydrofluoric acid																																	
Lacquer																																	
Liqueur etc.												✓	✓																				
Long length >150 mm		✓	✓																	✓													✓
Low ionic strength												✓	✓																				
Microsamples																												✓					✓
Milk	✓												✓												✓		✓						
Non-aqueous media																																	
Oil																																	
Paint													✓																				
Paper										✓																							
Penetration (cheese)																						✓					✓						
Photo baths												✓																					
Sea water, swimming pools	✓	✓	✓	✓			✓	✓	✓				✓		✓	✓	✓						✓	✓									
Soil							✓																										
Solids, high content of													✓																				
Suspensions												✓																					
Temperature, high or varying	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓																		
Temperature sensor, built-in								✓															✓	✓			✓						
TRIS buffer												✓	✓																				
Viscosity, high													✓																				
Voltammetry																																	
Yogurt, curdled milk							✓					✓	✓												✓		✓						

✓ Recommended     
 ■ Acceptable     
 ■ Not recommended     
 ■ Risk of damage

✓ There is a Radiometer Analytical electrode for every application. The above table gives examples of the many applications for which Radiometer Analytical electrodes can be used. Where possible, choose an electrode which is "Recommended" but where one electrode is required for several applications, it may be necessary to select one which is "Acceptable".

✓ In addition to your application, there are other considerations when choosing an electrode. The length and diameter obviously depend on the sample size and sample vessel. Radiometer Analytical offers electrodes with various lengths for specialist applications and small diameters for microsamples.

Full specifications of the electrodes featured in the table are given on the following pages.



Applications	Glass electrodes						Reference electrodes																		
	pHG200	pHG201	pHG211	pHG301	pHG311	XG250	REF200	REF201	REF251	REF321	REF361	REF401	REF421	REF451	REF601	REF621	REF921	XR110	XR150	XR200	XR400	XR440	XR820	MC20955b	
Alkaline solutions			✓		✓	✓		✓	✓	✓													✓	✓	
Aqueous solutions	✓	✓	✓	✓			✓	✓	✓	✓		✓	✓												
Blood		✓		✓																					
Education purposes				✓									✓						✓						
Emulsions		✓		✓																					
Fat/cream/cosmetics		✓		✓																					
Field use																									✓
Gel electrode																									
Hops/beer		✓		✓				✓		✓															
Hydrofluoric acid																									✓
Lacquer		✓		✓											✓										
Liqueur etc.		✓		✓																					✓
Long length >150 mm																									
Low ionic strength			✓		✓																				
Microsamples	✓						✓																		
Milk		✓		✓				✓		✓															
Non-aqueous media		✓		✓	✓						✓							✓							
Oil		✓		✓											✓										
Paint		✓		✓											✓										
Paper																									
Penetration (cheese)																									
Photo baths		✓		✓					✓						✓	✓	✓					✓			
Sea water, swimming pools	✓	✓		✓				✓		✓															
Soil																									
Solids, high content of		✓		✓																					✓
Suspensions		✓		✓																					
Temperature, high or varying	✓	✓	✓					✓	✓	✓															✓
Temperature sensor, built-in																									
TRIS buffer		✓		✓																					
Viscosity, high						✓																			
Voltammetry																					✓	✓	✓	✓	✓
Yogurt, curdled milk		✓		✓																					

✓ When measuring at high temperatures, the electrochemical system within the electrode needs to be considered. Radiometer Analytical offers Ag/AgCl electrodes specially designed for high-temperature measurements. See pages 6 and 7.

✓ To avoid chloride interference, a reference electrode with a chloride-free reference system must be used (e.g. Hg/Hg<sub>2</sub>SO<sub>4</sub> with K<sub>2</sub>SO<sub>4</sub> salt-bridge) or a reference electrode with a double salt-bridge construction. See pages 8 and 9.

## Combined pH Electrodes


To meet your analysis needs, combined pH electrodes are available with a variety of lengths, diameters and reference systems. To discover the advantages of Red Rod reference systems, please refer to page 3.

All our electrodes are delivered with a Certificate of Conformity signed by the Total Quality Manager. This specifies the sensitivity, zero pH, response time and batch number of the electrode.

We supply a bottle of filling solution with every combined pH electrode (with the exception of the GK2401B and GK2401C, and gel-filled electrodes).


You may need a cable or plug adapter to connect your chosen electrode to your meter. Please see page 19.

## Red Rod combined pH electrodes

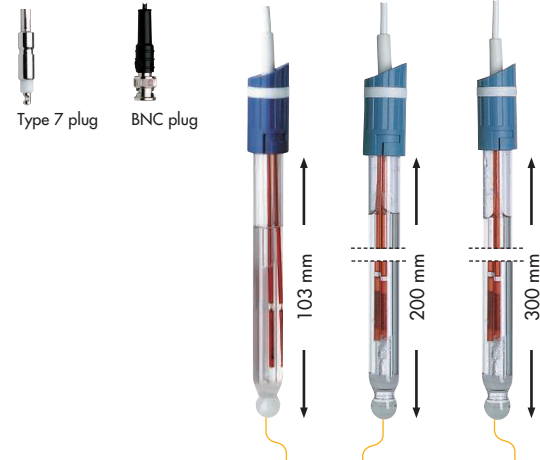

Applications	General		
Type	pHC2001	pHC2002	pHC2003
Part no. Type 7 -7			
Part no. BNC -8	E16M313	E16M315	E16M334
Reference system	Red Rod	Red Rod	Red Rod
pH range	0 - 12	0 - 12	0 - 12
Temperature range	-10 - 100°C	-10 - 100°C	-10 - 100°C
Diameter	12 mm	12 mm	12 mm
Minimum sample depth	18 mm	18 mm	18 mm
Liquid junction	Porous pin	Porous pin	Porous pin
Special features		Long	Long
Salt-bridge solution	Sat. KCl	Sat. KCl	Sat. KCl

## Combined pH electrodes

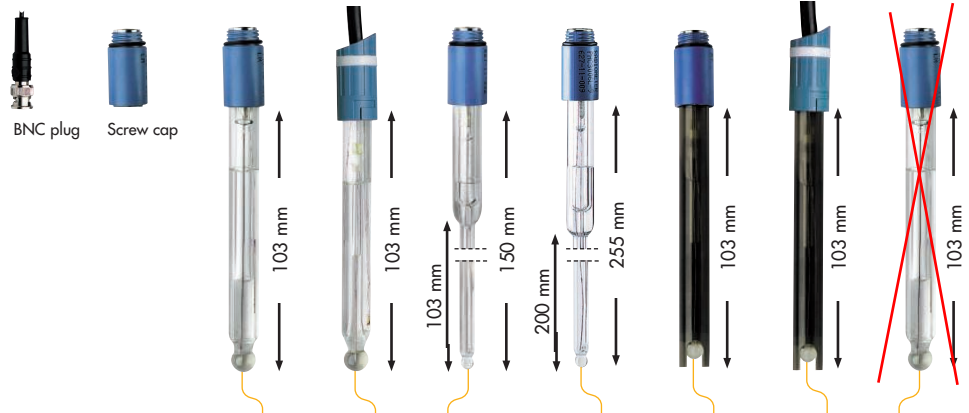

Applications	General-purpose						Alkaline samples
Type	pHC3001	pHC3081	pHC3006	pHC3006L	pHC3005	pHC3085	pHC3011
Part no. BNC -8	E16M300	E16M305			E16M302	E16M342	
Part no. Screw cap -9	E16M301		E16M306	E16M307	E16M303		E16M324
Reference system	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl
pH range	0 - 12	0 - 12	0 - 12	0 - 12	0 - 12	0 - 12	0 - 14
Temperature range	0 - 80°C	0 - 80°C	0 - 80°C	0 - 80°C	0 - 80°C	-5 - 80°C	0 - 80°C
Lower diameter	12 mm	12 mm	6.5 mm	6.5 mm	12 mm	12 mm	12 mm
Minimum sample depth	18 mm	18 mm	14 mm	14 mm	12 mm	12 mm	18 mm
Liquid junction	Porous pin	Porous pin	Porous pin	Porous pin	Porous pin	Porous pin	Porous pin
Special features		Temp. sensor	Small Ø	Long, small Ø	Robust	Temp. sensor	
Salt-bridge solution	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl

## Obtaining accurate and reproducible results: some tips

✓ Remember to make regular calibrations to ascertain the actual sensitivity of the electrode. Ideal electrode sensitivity lies within the range 97 - 100%. However, the electrode may be used within as wide a range as 95 - 10%.

✓ The calibration buffers should have the same temperature as the sample. Under normal conditions, a deviation of  $\pm 5^\circ\text{C}$  between sample and calibration buffers is acceptable.

✓ It is important to clean your electrodes regularly to ensure optimal response time. The GK•ANNEX Electrode Maintenance Kit is ideal for electrodes with a saturated KCl salt-bridge. See page 1.



purpose			Alkaline samples		Surface measurement	Soils	Cl <sup>-</sup> sensitive samples	Viscous samples	High KCl flow	Alkaline samples	General-purpose
pHC2401	pHC2005	pHC2085	pHC2011	pHC2015	pHC2441	pHC2051	pHC2501	pHC2601	pHC2701	GK2401B	GK2401C
E16M336	<del>E16M337</del>								<del>E16M339</del>	945-261	945-252
E16M400	E16M500	E16M501	E16M317	E16M318	E16M320	<del>E16M319</del>	<del>E16M321</del>	E16M322	E16M323		
Red Rod	Red Rod	Red Rod	Red Rod	Red Rod	Red Rod	<del>Red Rod</del>	<del>Red Rod</del>	Red Rod	Red Rod	Red Rod	Red Rod
0 - 12	0 - 12	0 - 12	0 - 14	0 - 14	0 - 12	<del>0 - 12</del>	<del>0 - 12</del>	0 - 12	0 - 12	0 - 14	0 - 12
-10 - 100°C	-10 - 100°C	-10 - 100°C	0 - 100°C	0 - 100°C	-10 - 100°C	<del>-10 - 100°C</del>	<del>-10 - 100°C</del>	-10 - 100°C	-10 - 100°C	0 - 100°C	-10 - 100°C
12 mm	12 mm	12 mm	12 mm	12 mm	12 mm	<del>12 mm</del>	<del>12 mm</del>	12 mm	12 mm	9.5 mm	9.5 mm
14 mm	14 mm	14 mm	14 mm	14 mm	14 mm	<del>18 mm</del>	<del>18 mm</del>	18 mm	18 mm	16 mm	16 mm
Annular ring	Porous pin	Porous pin	Porous pin	Porous pin	Annular ring	<del>Porous pin</del>	<del>Double/Porous</del>	Sleeve	Annular ring	Porous pin	Porous pin
Sat. KCl	Robust	Temp. sensor	Sat. KCl	Robust	Sat. KCl	<del>Sat. KCl</del>	<del>Sat. KCl</del>	Sat. KCl	Sat. KCl	<del>Sat. KCl</del>	Sat. KCl

Penetration	Microsamples		General-purpose Gel-filled		Penetration Gel-filled	General-purpose glass tubes		
pHC3031	XC161	pHC3359	pHC3105	pHC3185	pHC3131	pHC4000	pHC4001*	pHC4006*
E16M325	B10C161	E16M343	E16M308	E16M329	<del>E16M328</del>	<del>E16M331</del>	E16M332	<del>E16M323</del>
Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	<del>Ag/AgCl</del>	<del>Colomel</del>	Colomel	<del>Colomel</del>
0 - 12	0 - 12	0 - 12	2 - 12	2 - 12	<del>2 - 12</del>	<del>0 - 12</del>	0 - 12	<del>0 - 12</del>
0 - 80°C	0 - 80°C	0 - 80°C	0 - 80°C	0 - 80°C	<del>0 - 80°C</del>	<del>0 - 80°C</del>	0 - 80°C	<del>0 - 80°C</del>
8 mm	5 mm	3 mm	12 mm	12 mm	<del>8 mm</del>	<del>5 mm</del>	12 mm	<del>6.5 mm</del>
14 mm	8 mm	2 mm	18 mm	18 mm	<del>14 mm</del>	<del>9 mm</del>	8 mm	<del>4 mm</del>
Porous pin	Porous pin	Fibre	Opening	Opening	<del>Opening</del>	<del>Porous pin</del>	Porous pin	<del>Porous pin</del>
3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	3 M KCl + sat. AgCl	Solid gel (3 M KCl)	Temp. sensor	<del>Solid gel (3 M KCl)</del>	<del>Small Ø</del>	Sat. KCl	<del>Long, small Ø</del>

✓ Never touch the electrode bulb with your fingers. Any grease may affect the electrode membrane and cause a drifting potential.

✓ The typical response time for a combined pH electrode is 20 s depending on the sample, sample temperature, stirring etc.

✓ The expected lifetime of a combined pH electrode is 1<sup>1/2</sup> years provided that it has been correctly used and maintained.

\* These Electrodes Contain Mercury. Dispose According To Local, State Or Federal Regulations.

## Reference Electrodes


A reference electrode provides a stable potential whatever the measurement conditions. The main differences between reference electrodes are the type of reference system and the liquid junction.

Radiometer Analytical offers electrodes using both Red Rod and traditional technology (calomel, silver-silver chloride, mercurous sulphate and mercuric oxide).

Electrodes with a double junction are mostly designed for measurements using ion-selective electrodes or samples sensitive to  $\text{Cl}^-$  or  $\text{K}^+$ .


Sleeve junctions have a high flow to prevent the risk of clogging. Porous pins are used for general applications. Fibre reference electrodes have a slow flow to avoid contamination of the media.

## Red Rod reference electrodes

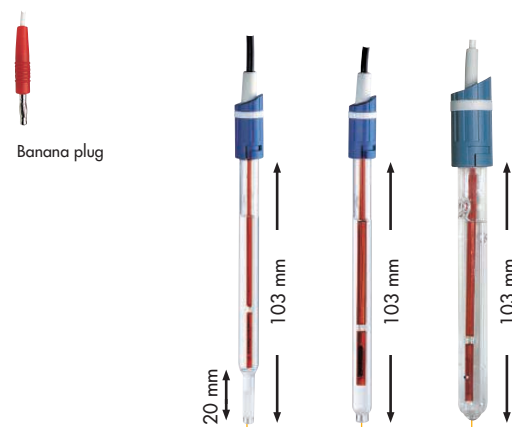

Applications	Micro-samples	General-purpose	
Reference system	Red Rod	Red Rod	Red Rod
Type	REF200	REF201	REF251
Part no.	E21M008	E21M009	E21M001
Temperature range	-10 – 100°C	-10 – 100°C	-10 – 100°C
Lower stem diameter	4 mm	7.5 mm	12 mm
Liquid junction	Porous pin	Porous pin	Porous pin/ porous pin
Special features	For use with pHG200		Double junction
Salt-bridge solution	Sat. KCl	Sat. KCl	Sat. KCl

## Reference electrodes

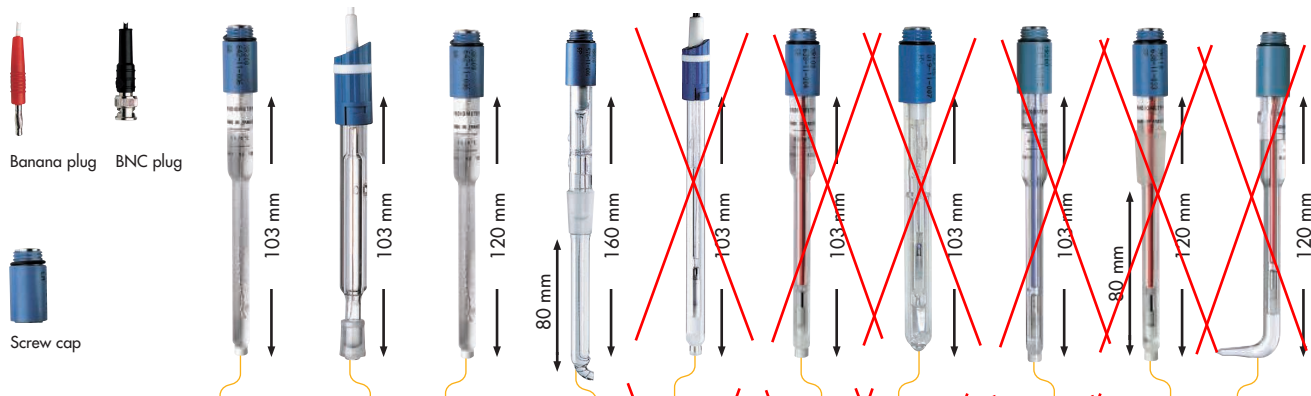


Applications	General-purpose	Non-aqueous media	Voltammetry	Impedance	Special-purpose			Non-aqueous media	Voltammetry
Reference system	Ag/AgCl	Ag/AgCl	Ag/AgCl	Ag/AgCl	Calomel	Calomel	Calomel	Calomel	Calomel
Type	REF321	REF361	XR300	XR820	REF401*	REF421*	REF451*	REF921*	XR110*
Part no.	E21M002	E21M003	B20B300	B20C820	E21M011	E21M004	E21M005	E21M007	B20B170
Temperature range	-5 – 80°C	-5 – 80°C	-5 – 80°C	-5 – 80°C	-10 – 60°C	-10 – 60°C	-10 – 60°C	0 – 80°C	-10 – 60°C
Lower diameter	8 mm	12 mm	8 mm	8 mm	7.5 mm	8 mm	12 mm	8 mm	8 mm
Liquid junction	Porous pin	Glass sleeve	Porous pin	Porous pin/ fibre	Porous pin	Porous pin	Porous pin/ porous pin	Porous pin	Porous pin
Connection	Screw cap	BNC plug shielded cable	Screw cap	Screw cap	Banana plug	Screw cap	Screw cap	Screw cap	Screw cap
Salt-bridge solution	3 M KCl with sat. AgCl	3 M KCl with sat. AgCl	Sat. KCl + sat. AgCl	Sat. KCl + sat. AgCl upper res	Sat. KCl	Sat. KCl	Sat. KCl in inner reservoir	1 M LiCl	Sat. KCl
									Sat. KCl






## Salt-bridge junctions

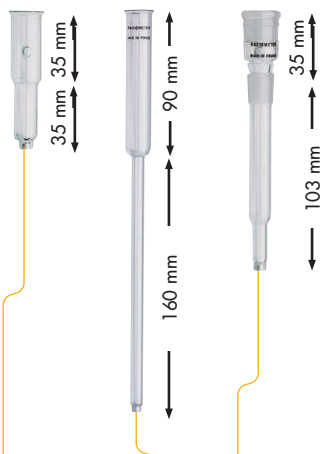

Type	AL100	AL110	AL120
Part no.	B40A520	B40A550	B40A610
Lower diameter	12 mm	8 mm	8 mm
Junction	Ceramic	Ceramic	Ceramic
Special features	for Ø 7.5 or 8 mm ref. electrodes		for ref. electrodes with ground joint

During measurement, the salt-bridge junction (also called liquid junction protection tube) prevents the reference electrode from being contaminated by the media and vice-versa. It consists of a glass tube terminated by a liquid junction. Always ensure the filling solution used corresponds to the application.

## Temperature sensors

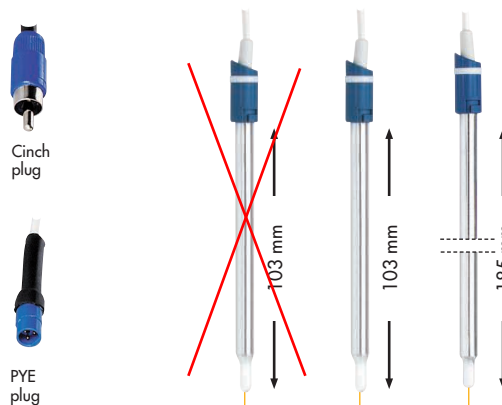

Type	<del>T101</del>	T201	T202
Part no.	<del>E51M003</del>	E51M001	E51M002
Temperature range	<del>-10 – 105°C</del>	-10 – 105°C	-10 – 105°C
Diameter	<del>7.5 mm</del>	7.5 mm	7.5 mm
Material	<del>Glass</del>	Glass	Glass
Connection	<del>PYE</del>	Cinch	Cinch
Instrument	<del>PHM83/85 ON83/85 CDM83</del>	PHM2XX CDM2XX PHM9X	PHM2XX CDM2XX PHM9X


For pH measurements and measurements with ion-selective electrodes and conductivity cells, use of a temperature sensor is recommended because temperature affects the results.

## Obtaining accurate and reproducible results: some tips

✓ To ensure a stable reference potential, the salt-bridge solution should not be more than 5 mm below the electrode filling hole. Top up when required.

✓ Always read the proper storage conditions described in the operating instructions.

✓ Clean your electrodes regularly. The operating instructions give details of the most suitable cleaning procedure. A correctly maintained reference electrode has a lifetime of approximately 2 years.

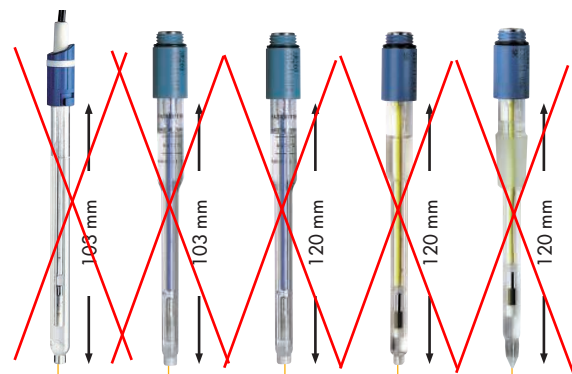

Chloride free		Voltammetry		
Hg/Hg <sub>2</sub> SO <sub>4</sub>	Hg/Hg <sub>2</sub> SO <sub>4</sub>	Hg/Hg <sub>2</sub> SO <sub>4</sub>	Hg/HgO	Hg/HgO
<del>REF601*</del>	<del>REF621*</del>	<del>XR200*</del>	<del>XR400*</del>	<del>XR440*</del>
<del>E21M012</del>	<del>E21M006</del>	<del>B20B200</del>	<del>B20B400</del>	<del>B20B440</del>
<del>0 – 60°C</del>	<del>0 – 60°C</del>	<del>0 – 60°C</del>	<del>0 – 60°C</del>	<del>0 – 60°C</del>
<del>7.5 mm</del>	<del>8 mm</del>	<del>8 mm</del>	<del>8 mm</del>	<del>8 mm</del>
<del>Porous pin</del>	<del>Porous pin</del>	<del>Porous pin</del>	<del>Porous pin</del>	<del>Fibre</del>
<del>Banana</del>	<del>Screw cap plug</del>	<del>Screw cap</del>	<del>Screw cap</del>	<del>Screw cap</del>
<del>Sat. K<sub>2</sub>SO<sub>4</sub></del>	<del>Sat. K<sub>2</sub>SO<sub>4</sub></del>	<del>Sat. K<sub>2</sub>SO<sub>4</sub></del>	<del>0.1 M KOH</del>	<del>1 M KOH</del>

\* These Electrodes Contain Mercury.  
Dispose According To Local, State Or Federal Regulations.

## Glass pH Electrodes


pH glass electrodes from Radiometer Analytical are made with highly shock-resistant glass, making them extremely tough. The average lifetime of a correctly maintained glass pH electrode is approximately 2 years.

Certain electrodes are made with special glass for pH measurements in strongly alkaline solutions. These electrodes allow measurements within the 0 - 14 pH range without significant deviation from the theoretical response.

Glass pH electrodes are always used in combination with reference electrodes. To find the right reference electrode for your application, refer to pages 8 and 9.


The zero pH of Red Rod glass pH electrodes when used with their recommended reference electrode is approximately pH 6.65. The zero pH of glass pH electrodes is approximately pH 6.65 when measured against a calomel reference electrode (e.g. REF421, XR110) and approximately pH 7.25 when measured against an Ag/AgCl reference electrode (e.g. REF321).

### Red Rod glass pH electrodes

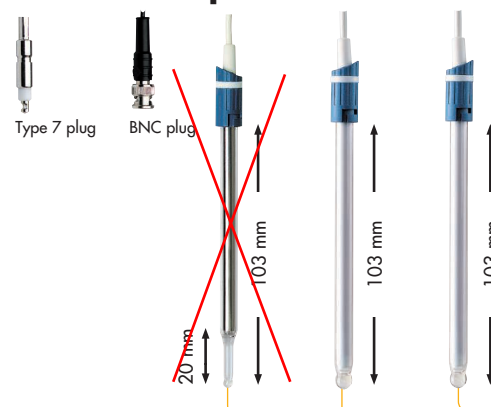

Applications	<del>Micro-samples</del>	General-purpose	Alkaline samples
Type	<del>pHG200</del>	pHG201	pHG211
Part no. Type 7		E11M005	
Part no. BNC	E11M002	E11M006	E11M003
pH range	0 - 12	0 - 12	0 - 14
Temperature range	10 - 100°C	-10 - 100°C	10 - 100°C
Lower diameter	4 mm	7.5 mm	7.5 mm
Reference system	Red Rod	Red Rod	Red Rod
Minimum sample depth	4 mm	7 mm	7 mm
Recommended reference electrode	REF200	REF201/251/261 REF201 for general use	

## Glass pH electrodes

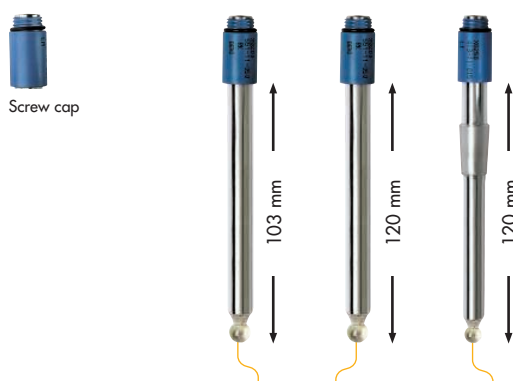


Applications	General-purpose	Alkaline samples	
Type	pHG301	pHG311	XG250
Part no.	E11M001	E11M004	B10B250
pH range	0 - 12	0 - 14	0 - 14
Temperature range	0 - 80°C	0 - 80°C	0 - 80°C
Diameter	12 mm	12 mm	10 mm*
Reference system	Ag/AgCl	Ag/AgCl	Ag/AgCl
Minimum sample depth	12 mm	12 mm	12 mm
Recommended reference electrodes	REF421 REF321	REF421 REF321	XR100


\* under ground joint

## Metal and Combined Metal Electrodes

Single and combined metal electrodes can be used for redox measurements, redox titrations or potentiometric techniques whereas double metal electrodes are ideal for imposed current potentiometric titration. Radiometer Analytical's electrode range features single, double and combined metal electrodes with platinum, silver, gold, antimony, mercury and glassy carbon sensing elements.

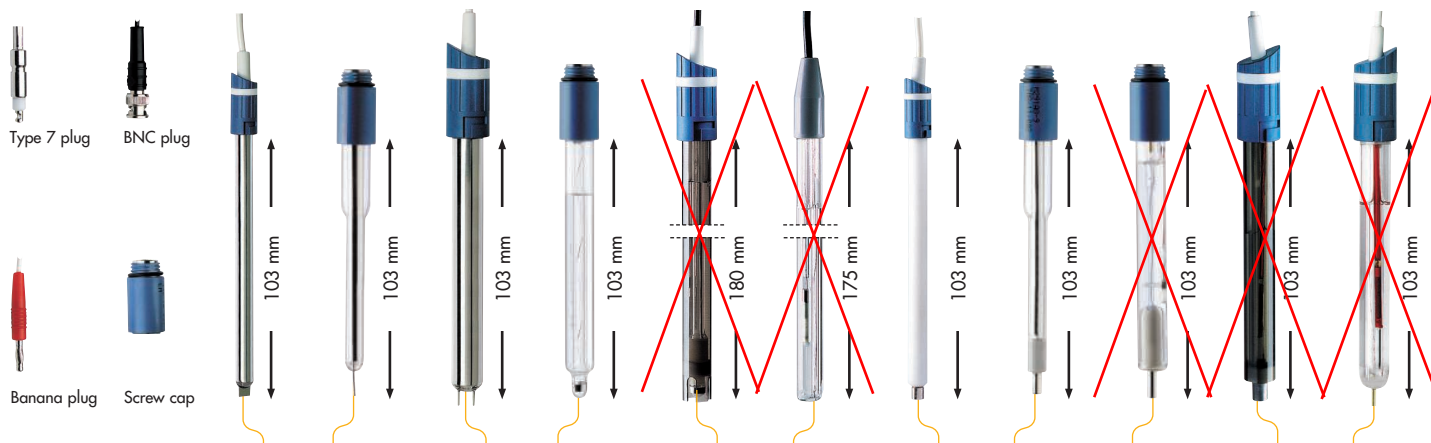

Applications	General-purpose		Karl Fischer titrations	Redox meas.	Kappa no. titrations	COD titrations	Silver halides Titration			pH meas. titration	Redox meas.
Type	M241Pt	M231Pt	M231Pt2	MC3051Pt	MC408Pt*	MC602Pt*	M295Ag	M291Ag	MC6091Ag	MC2095Sb	MC201Au-8
Part no.	E31M001	E31M002	E32M001	E31M003	<del>E45-390</del>	<del>E45-360</del>	E34M003	E34M002	<del>E34M004</del>	E36M001	<del>E36M005</del>
Sensing element	Platinum plate 5 x 5 mm	Platinum wire Ø 1 mm	2 platinum wires Ø 1 mm	Platinum ring	<del>Platinum ring</del>	<del>Platinum wire</del>	Silver rod	Silver rod	<del>Silver rod</del>	Antimony rod	<del>Gold wire Ø 1 mm</del>
Temperature range	-10 - 100°C	-10 - 100°C	-10 - 100°C	0 - 80°C	<del>-10 - 60°C</del>	<del>0 - 60°C</del>	-10 - 100°C	0 - 80°C	<del>0 - 60°C</del>	-10 - 100°C	<del>-10 - 100°C</del>
Lower diameter	7.5 mm	8 mm	12 mm	12 mm	<del>12 mm</del>	<del>9.5 mm</del>	7.5 mm	8 mm	<del>12 mm</del>	12 mm	<del>12 mm</del>
Reference system				Ag/AgCl	<del>Hg/Hg<sub>2</sub>Cl<sub>2</sub></del>	<del>Hg/Hg<sub>2</sub>SO<sub>4</sub></del>			<del>Hg/Hg<sub>2</sub>SO<sub>4</sub></del>	Red Rod	<del>Red Rod</del>
Liquid junction				Porous pin	<del>Porous pin</del>	<del>Porous pin</del>			<del>Porous pin</del>	Red Rod	<del>Red Rod</del>
Connection	Banana plug	Screw cap	2 x banana	Screw cap	Type 7 plug	Type 7 plug	Banana plug	Screw cap	<del>Screw cap</del>	BNC plug	<del>BNC plug</del>
Salt-bridge solution				3 M KCl with sat. AgCl	<del>Sat. KCl</del>	<del>Sat. K<sub>2</sub>SO<sub>4</sub></del>			<del>Sat. K<sub>2</sub>SO<sub>4</sub></del>	Sat. KCl	<del>Sat. KCl</del>

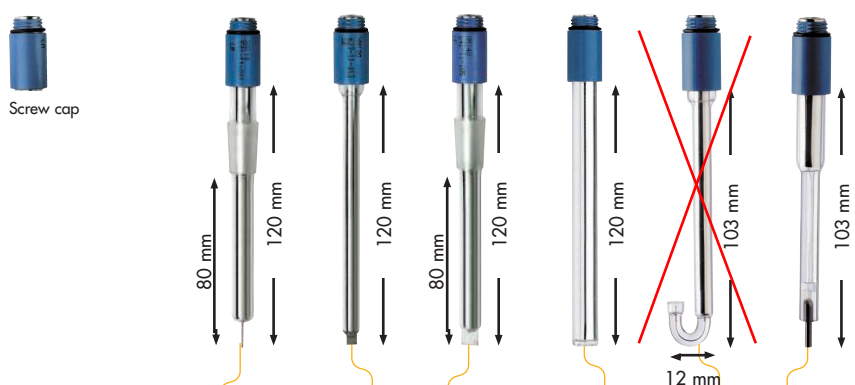

Applications	Voltammetry				Electro-chemistry	Voltammetry Redox
Type	XM110	XM120	XM140	XM150	M221Hg*	M291C
Part no.	B35M110	B35M120	B35M140	B35M150	<del>E36M002</del>	E36M003
Sensing element	Platinum wire Ø 1 mm	Platinum plate 5 x 5 mm	Platinum plate 8 x 8 mm	Platinum disc Ø 10 mm	<del>Mercury cup</del>	Glassy carbon rod Ø 3 mm
Temperature range	-10 - 100°C	-10 - 100°C	-10 - 100°C	-10 - 100°C	<del>-10 - 100°C</del>	0 - 60°C
Diameter	8 mm	8 mm	8 mm	12 mm	<del>8 mm</del>	8 mm
Connection	Screw cap	Screw cap	Screw cap	Screw cap	<del>Screw cap</del>	Screw cap

\* These Electrodes Contain Mercury. Dispose According To Local, State Or Federal Regulations.

\*\* This Electrode Requires Mercury To Function. Dispose According To Local, State Or Federal Regulations.


## Ion-selective Electrodes

The use of ion-selective electrodes constitutes a quick and cost-effective method of determining ionic activity.

Due to their quality and diversity, Radiometer Analytical ion-selective electrodes efficiently solve various analytical problems of species in solution.

For details of the reference electrodes recommended for use with ion-selective electrodes, please see pages 8 and 9.

Each electrode is carefully checked before packing and is delivered with an individual control report including information such as the response curve, checkpoint value and slope.

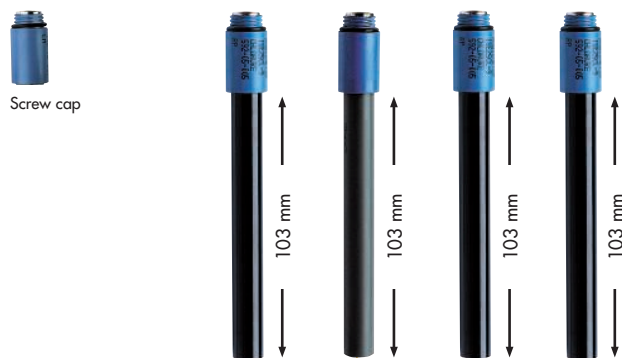


Ion Species	Br <sup>-</sup>	Ca <sup>++</sup>	Cl <sup>-</sup>	Cl <sup>-</sup>
Type	ISE25Br	ISE25Ca	ISE25Cl	ISE/HS25Cl
Part no.	E41M001	E41M002	E41M003	E41M004
Concentration range (M)	10 <sup>-6</sup> – 10 <sup>0</sup>	2 x 10 <sup>-6</sup> – 10 <sup>0</sup>	5 x 10 <sup>-5</sup> – 10 <sup>0</sup>	10 <sup>-6</sup> – 10 <sup>0</sup>
Concentration range (ppm)	0.08 – 80000	0.1 – 40000	1.8 – 35000	0.04 – 35000
Temperature range	0 – 60°C	0 – 50°C	0 – 60°C	0 – 60°C
Diameter	12 mm	12 mm	12 mm	12 mm
pH range	2 – 14	4 – 12	2 – 14	2 – 4
Sensor type	Solid-state	PVC membrane	Solid-state	Solid-state
Interfering ions	I <sup>-</sup> , S <sup>2-</sup> , CN <sup>-</sup> , SCN <sup>-</sup>	Zn <sup>++</sup> , Pb <sup>++</sup>	I <sup>-</sup> , CN <sup>-</sup> , Br <sup>-</sup> , SCN <sup>-</sup>	I <sup>-</sup> , CN <sup>-</sup> , Br <sup>-</sup> , SCN <sup>-</sup>
Recommended reference electrodes	REF251	REF201/REF251 REF451	REF251 REF451 REF601	REF251 REF451 REF601
Inner solution		S41M001		
Membrane kit part no. [2]		E91M001		

\* Equally sensitive to Ag<sup>+</sup>.

(1) Any organic anions or cations that chemically resemble the species being analysed.  
 (2) The membrane kit contains three electrode tubes and a bottle of inner solution. The - M27Ag-9 (part no. E34M001). Electrode insert with screw cap.

Note: REF251 and REF451 are double-junction electrodes. You may need to change the

### Analysis with ion-selective electrodes is reliable and convenient when a few simple guidelines are followed:

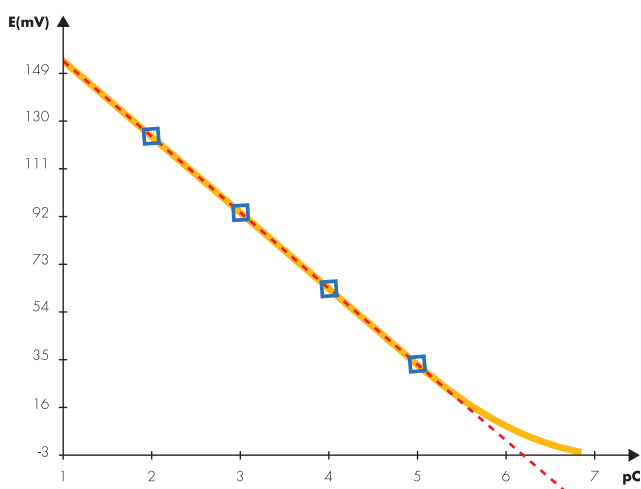

Typical calibration curve for an ion-selective electrode.

✓ Ion-selective electrodes respond logarithmically over several decades of concentration. However, to be able to rely on the result in the area near the detection limit, you must use several calibration points - see figure.

✓ Low calibration solutions should be prepared each time the electrode is calibrated.

✓ Always calibrate the lowest standard concentration first.

✓ Make sure that calibration standards and samples have the same pH value, temperature and ionic strength. If necessary, add ionic strength adjustment buffer (TISAB/ISA solutions).

✓ Some ion-selective electrodes are sensitive to light so use black beakers (part no. 904-515) when measuring Cl<sup>-</sup>, I<sup>-</sup> and Br<sup>-</sup>.




103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm	103 mm
<b>CN<sup>-</sup></b>	<b>Cu<sup>++</sup></b>	<b>F<sup>-</sup></b>	<b>F<sup>-</sup></b>	<b>I<sup>-</sup></b>	<b>K<sup>+</sup></b>	<b>Na<sup>+</sup></b>	<b>NH<sub>4</sub><sup>+</sup></b>	<b>NO<sub>3</sub><sup>-</sup></b>	<b>Pb<sup>++</sup></b>	<b>S<sup>-</sup>(Ag<sup>+</sup>)<sup>+</sup></b>	<b>Surfactant</b>
<b>ISE25CN</b>	<b>ISE25Cu</b>	<b>ISE25F</b>	<b>ISEC301F</b>	<b>ISE25I</b>	<b>ISE25K</b>	<b>ISE21Na</b>	<b>ISE25NH4</b>	<b>ISE25NO3</b>	<b>ISE25Pb</b>	<b>ISE25S</b>	<b>ISE25X</b>
E41M005	E41M006	E41M007	E41M017	E41M008	E41M009	E41M010	E41M013	E41M014	E41M015	E41M016	E41T001
$5 \times 10^{-7} - 10^{-3}$	$10^{-6} - 10^0$	$5 \times 10^{-7} - 10^0$	$5 \times 10^{-7} - 10^0$	$10^{-6} - 10^{-1}$	$2 \times 10^{-6} - 10^0$	$2 \times 10^{-6} - 10^0$	$3 \times 10^{-6} - 10^0$	$3 \times 10^{-6} - 10^0$	$10^{-6} - 10^0$	$5 \times 10^{-7} - 10^0$	$5 \times 10^{-2} - 10^{-5}$
0.013 - 25	0.06 - 60000	0.01 - 20000	0.01 - 20000	0.13 - 13000	0.08 - 40000	0.04 - 20000	0.06 - 20000	0.2 - 60000	0.2 - 200000	0.02 - 30000	
0 - 60°C	0 - 60°C	0 - 60°C	0 - 60°C	0 - 60°C	0 - 50°C	0 - 70°C	0 - 50°C	0 - 50°C	0 - 60°C	0 - 60°C	0 - 50°C
12 mm	12 mm	7 mm	12 mm	12 mm	2 mm	12 mm	12 mm	12 mm	12 mm	12 mm	12 mm
9 - 14	0 - 14	5-11 (10 <sup>-1</sup> M F)	5-7 (10 <sup>-6</sup> M F)	2 - 12	2 - 12	5 - 12	3 - 8	3 - 10	3 - 7	11 - 14	2 - 12
		(3-7 direct meas.)	5-7 (10 <sup>-6</sup> M F)					(direct meas.)	(direct meas.)		
Solid-state	Solid-state	Solid-state	Solid-state	Solid-state	PVC membrane	Glass	PVC membrane	PVC membrane	Solid-state	Solid-state	PVC
I <sup>-</sup> , Br <sup>-</sup> , S <sup>2-</sup>	Hg <sub>2</sub> , Ag <sup>+</sup> , Cl <sup>-</sup>	OH <sup>-</sup>	OH <sup>-</sup>	CN <sup>-</sup> , Hg <sup>++</sup> , Br <sup>-</sup> , SCN <sup>-</sup>	NH <sub>4</sub> <sup>+</sup> , Rb <sup>+</sup>	Li <sup>+</sup> , NH <sub>4</sub> <sup>+</sup> , Rb <sup>+</sup> , Ag <sup>+</sup>	K <sup>+</sup>	Cl <sup>-</sup> , Br <sup>-</sup> , NO <sub>2</sub> <sup>-</sup>	Hg <sup>++</sup> , Cu <sup>++</sup> , Ag <sup>+</sup>	Precipitated cations forming complexes with sulphides	See note (1)
REF201/REF251 REF451	REF201/REF251 REF451	REF201/REF251 REF451	Built-in Ag/AgCl reference	REF201/REF251 REF451	REF251 REF451	REF251 REF451	REF251 REF451	REF201/REF251 REF451	REF251 REF451	REF251 REF451	REF251
					S41M002 E91M002		S41M003 E91M003	S41M004 E91M004			

electrode insert can be ordered separately:

filling solution in the second compartment for certain measurements (see operating instructions).


## Electrode membranes

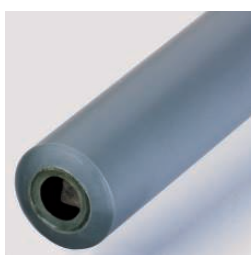

Calcium



Sodium

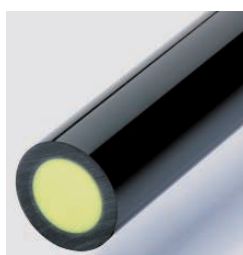

Cyanide



Fluoride

The sensing element of an ion-selective electrode is the membrane. The type of membrane varies depending on the species to be measured. Membranes may be solid-state (e.g. cyanide or fluoride), PVC (e.g. calcium) or glass (sodium).

### Solutions for fluoride ion selective electrodes

Ion and Electrode	Description	Part no.
F <sup>-</sup> ISE25F	F 1000 : 1000 ppm standard solution, Fluoride, 500 ml	S41M007
F <sup>-</sup> ISEC301F	RENOVO.F, cleaning solution for ISE electrodes, 500 ml	S41M008
F <sup>-</sup> ISEC25F	<del>K3120-30, KCl &amp; AgCl saturated solution, 30 ml</del>	<del>S21M009</del>
F <sup>-</sup> ISEC301F		

**S21M004**  
**S21M011**

## Conductivity Cells

It is important to choose a cell with the right construction and geometry for your particular application and working conditions. Radiometer Analytical offers conductivity cells for a wide variety of applications.

2-pole cells have a traditional design based on two plates of platinum. They are ideal for routine measurement of conductivity and for use with a sample changer due to the easy rinsing.

3-pole cells consist of 3 platinum rings which facilitate optimal shielding during measurement.

4-pole cells consist of 4 platinum rings. They ensure accurate results over several decades of conductivity with a single cell using just one calibration. They are particularly recommended when performing high conductivity measurements.

With a 4-pole cell, an alternating current is applied to the two outer rings and the voltage is measured on the 2 inner rings thereby avoiding errors due to polarisation effects and guaranteeing measurement accuracy. The CDC566T and CDC866T Conductivity Cells offer the advantages of this 4-pole design with a built-in temperature sensor. The tough epoxy body can be easily removed for rinsing which makes these cells ideal for measurements across a wide conductivity range even in harsh conditions.

### Some tips

✓ Conductivity measurement is temperature dependent (if the temperature increases so does the conductivity value). When performing conductivity measurements, it is advisable to choose a cell with a built-in temperature sensor or use a separate sensor.

✓ Remember to calibrate your conductivity cell regularly as the cell constant may vary due to changes in electrode surface due to contamination, for example.

✓ If your cell is to be used with one of the conductivity meters in our current range, the user-friendly CDM210 or the high-performance CDM230, order a cell with a MAB6 plug. This plug type is also suitable for the CDM92 Conductivity Meter. For older meters (CDM80 or CDM83), choose a cell with a UHF plug.



Applications	General-purpose				
Type	CDC566T	CDC866T	CDC565	CDC641T	CDC745-9
Part no. UHF					
Part no. MAB6	E61M010	E61M015	E61M003	B15B001	
Part no. screw cap					E61M013
Body	Epoxy <sup>1)</sup>	Epoxy <sup>1)</sup>	Epoxy	Glass	Epoxy
Cell. constant (cm <sup>-1</sup> ) <sup>9)</sup>	1.0	1.0	1.0	1.0	1.0
Number of poles	4	4	4	2	2
Platinised	NO	YES	NO	YES	YES
Temperature sensor	YES	YES	NO	YES	NO
Diameter	12 mm	12 mm	12 mm	12 mm	12 mm
Max. temperature	80°C	80°C	80°C	100°C	100°C
Min. immersion depth	35 mm	35 mm	30 mm	14 mm	14 mm
CDM210/CDM230	YES	YES	YES	YES	YES <sup>3)</sup>
CDM92	YES	YES	YES	YES	YES <sup>3)</sup>
CDM80	NO	NO	NO	NO	NO
CDM83	NO	NO	NO	YES <sup>2)</sup>	YES <sup>4)</sup>

1) Removable epoxy tube for easy rinsing (can be replaced by glass tube - see accessories)

2) Use adapter part no. A94P002

3) Use cable part no. A94L136

4) Use cable part no. A94L119

5) Polymethylpentene. Removable part for easy rinsing

6) Minimum sample volume

7) Diameter below the head

8) Use adapter part no. A94P001

9) The cell constant is determined experimentally for each cell and the value is given as a guideline only

## Recommended conductivity cells by application

Applications/Characteristics	CDC566T	CDC866T	CDC565	CDC641T	CDC745-9	XE100	CDC511T	CDC741T	CDC861T	CDC267-9	CDC241-9	CDC749	CDC104	CDC304
Wide conductivity range (general purpose)	✓	✓	✓					✓						
Various aqueous and non-aqueous media			✓		✓			✓				✓	✓	
Very strong acids and bases								✓						
Use with sample changer							✓							
Built-in temperature sensor	✓	✓	✓			✓	✓	✓						
Continuous measurements			✓		✓					✓				
Microsamples											✓			
Flow measurements						✓								
Titration										✓				
Salinity (high conductivity)		✓						✓						
Pure water							✓							
Meets requirements of USP 24-NF19			✓			✓								
Meets requirements of EP 2.2.38				✓										
Use in glass tubes											✓			
Plastic body	✓	✓	✓		✓	✓								
Viscous media										✓				
Highly resistant media									✓					
Field use	✓	✓	✓		✓									



	Pure water	Sample changer	Strong acids, bases	Resistant media	Viscous Titration	Micro-samples	Immersion	
	CDC511T	CDC741T	CDC861T	CDC267-9	CDC241-9	CDC749	CDC104	CDC304
	E61M009	E61M012	E61M016			E61M014	945-000	945-002
B60E100				E61M011	E61M008			
Glass	TPX <sup>§1</sup>	Glass	Glass	Glass	Glass	Glass	Glass	Glass
1.0	1.0	1.0	1.0	0.1	1.0	1.7	1.0	1.0
2	4	2	4	2	2	2	3	3
YES	NO	YES	YES	NO	NO	YES	YES	NO
NO	YES	YES	YES	NO	NO	NO	NO	NO
12 mm	12 mm	12 mm	12 mm	12 mm	12 mm	4 mm	12 mm <sup>§2</sup>	12 mm <sup>§2</sup>
100°C	80°C	100°C	100°C	100°C	100°C	100°C	100°C	100°C
10 mm	3 ml <sup>§1</sup>	10 mm	35 mm	26 mm	10 mm	4 mm	55 mm	55 mm
YES <sup>§1</sup>	YES	YES	YES	YES <sup>§1</sup>	YES <sup>§1</sup>	YES	YES <sup>§1</sup>	YES <sup>§1</sup>
YES <sup>§1</sup>	YES	YES	YES	YES <sup>§1</sup>	YES <sup>§1</sup>	YES	YES <sup>§1</sup>	YES <sup>§1</sup>
NO	NO	NO	NO	NO	YES <sup>§1</sup>	NO	YES	YES
YES <sup>§1</sup>	NO	YES <sup>§1</sup>	NO	YES <sup>§1</sup>	YES <sup>§1</sup>	YES <sup>§1</sup>	YES	YES

Part no.	Accessory
X31M013	Epoxy tube for CDC566T/CDC866T, diameter 12 mm
X51M002	Flow cell for 12 mm diameter sensors
X31M014	Circulation/pipette piece for CDC511T with set of accessories

# Certified Standards

## Accredited Calibration Laboratory

The use of standards manufactured by an accredited laboratory gives you complete confidence in the traceability chain and calculated uncertainties. Radiometer Analytical is part of the Hach Lange group whose Calibration Laboratory, located in Berlin, is accredited for the calibration of certified pH and conductivity standards by DKD, the German national accreditation authority. This accreditation (No. DKD-K-47901) is recognised in more than 35 countries. All standards are formulated in compliance with NIST, IUPAC and DIN19266 specifications<sup>(1)</sup>.

IUPAC Series **certified pH standards** are fully traceable to the Standard Hydrogen Electrode through Certified Reference Material produced by the Primary Laboratory of either NIST or DFM<sup>(2)</sup>.

Radiometer Analytical **certified conductivity standards** are fully traceable to SI units via Certified Reference Materials manufactured by NIST. They are prepared and calibrated according to the internationally recognised demal scale laid down by the International Organisation of Legal Metrology (OIML – Recommendation No. 56).

## Totally reliable packaging

Each standard is delivered in an airtight can and is supplied with a film-wrapped Certificate of Conformity and Traceability as well as a DKD Calibration Certificate specifying the exact measured value and calculated expanded uncertainty ( $k=2$ ) for each batch of pH and conductivity solution.

**A guaranteed value**  
the nominal value of the standard is guaranteed until first opening, even after several years' storage<sup>(3)</sup>.

**A useful reference**  
a pH or conductivity temperature dependence table is printed on the bottle.

**Safe storage**  
thanks to the airtight can.

**GLP - right from the start**  
the date of opening can be marked directly on the bottle.

**A Proven traceability**  
each standard comes with its own Certificate of Conformity and Traceability.

**Recommended shelf life**  
the observed typical shelf life after opening is printed on the bottle, ensuring timely use (2-3 months depending on the standards).



**GLP - every day**  
the number of uses can be marked directly on the bottle in line with Good Laboratory Practice.

**Accurate values**  
the nominal value and the tolerance of pH standards are given with a resolution of 3 significant decimals.

**B International accreditation**  
a DKD Calibration Certificate proves the traceability to national standards.

**Foolproof calibration**  
each bottle of pH standard is supplied with 25 colour-coded beakers.

<sup>(1)</sup> Except pH7.000 formulated by us according to NIST, IUPAC and DIN19266 specifications of pH6.865 and pH7.413.

<sup>(2)</sup> Dansk Fundamental Metrologi A/S is accredited for pH measurements by the Danish Accreditation and Metrology Fund (DANAK) (no. 255).

<sup>(3)</sup> 4-year guarantee for pH standards (2 years for pH 12.45). 2-year guarantee for conductivity standards (1 year for the 25  $\mu\text{S}/\text{cm}$ ). The guarantee is valid from the date of the Certificate of Conformity and Traceability and only applies to standards kept unopened in the sealed can in compliance with the storage recommendations specified in the certificate.

<sup>(4)</sup> ISO Guide 31 provides guidelines for certified reference materials.

### Glossary & useful links

NIST	National Institute of Standards and Technology <a href="http://www.nist.gov">www.nist.gov</a>
IUPAC	International Union of Pure and Applied Chemistry <a href="http://www.iupac.org">www.iupac.org</a>
DFM	Danish National Metrology Institute <a href="http://www.dfm.dtu.dk">www.dfm.dtu.dk</a>
DKD	Deutscher Kalibrierdienst <a href="http://www.dkd.eu">www.dkd.eu</a>
DANAK	Danish Accreditation and Metrology Fund <a href="http://www.danak.org">www.danak.org</a>
OIML	International Organisation of Legal Metrology <a href="http://www.oiml.org">www.oiml.org</a>





**A**

**Certificate of Conformity and Traceability**  
**Zertifikat über Konformität und Rückführbarkeit**  
**Certificat de Conformité et de Traçabilité**

**pH 10.012 ±0.010 (k=2) – 25°C**

Part No.:  
Bestellnr.: S11M007  
Code :

Batch No.:  
Chargenr.: C01311  
Lot n° :

**Traceability:**  
Traceable to IUPAC/NIST pH scale.  
Certified reference material: Radiometer Medical A/S ampoule of pH standard, batch **49.48**. Nominal value pH 10.012 at 25°C. The exact value of this standard was determined with an expanded uncertainty of ±0.0025 pH by Radiometer Medical A/S's accredited Primary Laboratory using a standard Hydrogen Electrode Apparatus. Samples from the batch are stored at Hach Lange GmbH for the warranty period of 4 years from the date of issue of the certificate.

**Rückführbarkeit:**  
Rückführbar auf die IUPAC/NIST pH-Skala. Zertifikates-Referenzmaterial: Radiometer Medical A/S Präzisionsampullenlösung, Charge **49.48**, pH-Nominalwert 10.012 bei 25°C. Der exakte Wert dieses Standards wurde mit einer erweiterten Gesamtsicherheits von ±0.0025 pH in dem akkreditierten Primärlaboratorium von Radiometer Medical A/S unter Verwendung einer Standard-Wasserstoff-Elektrodenvorrichtung bestimmt. Rückstell:

**Tracabilité :**  
Traçable selon l'échelle pH de IUPAC/NIST. Matériau de référence certifié : étalon pH Radiometer Medical A/S en ampoule, lot **49.48**. Valeur nominale pH 10.012 à 25°C. La valeur exacte de cet étalon a été déterminée avec une incertitude élargie de ±0.0025 pH par le Laboratoire Primaire accrédité de Radiometer Medical A/S à l'aide d'un système utilisant une électrode standard à hydrogène. Des échantillons issus du matériel, lot sont

**Certified value:**  
pH 10.012 ±0.010 (k=2) at 25°C  
The limits of the expanded uncertainty are given to guarantee a confidence level of ~95% (k=2). This uncertainty reflects the combined effects of measurement errors, variability among bottles and possible changes during storage in unopened tin.

**Stability:**  
When stored in an unopened tin, the certified value is guaranteed for 4 years from the date of issue of the certificate.

**Homogeneity:**  
7 bottles were selected for analytical control. Results from different bottles showed no statistically significant differences, nor was there any correlation between values obtained and the bottling sequence.

**Measurement:**  
The certified value was determined by measurements of samples with dedicated electrodes under thermostated conditions using a high-resolution meter (0.0004 pH) traceable to electrical primary standards. Separate bottles were controlled for bacterial and mould contamination before the batch was released.

**Date / Datum / Date :** 2007-05-23  
TTTTMMDD / JJJJMMJJ / AAAA.MM.JJ  
**Signature / Unterschrift / Signature :**

F014K03H

**Full and recognised traceability**

Each bottle of certified pH or conductivity standard manufactured by Radiometer Analytical is supplied with two certificates: a Certificate of Conformity and Traceability drawn up according to ISO Guide 31<sup>(4)</sup> and a DKD Calibration Certificate. The back of the first certificate contains information about the composition and preparation of the standard as well as recommendations for storage and use to ensure optimum performance.

**B**

**IUPAC Series certified pH standards**

Type	Value	Qty	Part No.
pH1.679	pH 1.679 ±0.010 at 25°C	500 ml	S11M001
pH4.005	pH 4.005 ±0.010 at 25°C	500 ml	S11M002
pH6.865	pH 6.865 ±0.010 at 25°C	500 ml	S11M003
pH7.000*	pH 7.000 ±0.010 at 25°C	500 ml	S11M004
pH7.413	pH 7.413 ±0.010 at 25°C	500 ml	S11M005
pH9.180	pH 9.180 ±0.010 at 25°C	500 ml	S11M006
pH10.012	pH 10.012 ±0.010 at 25°C	500 ml	S11M007
pH12.45	pH 12.45 ±0.05 at 25°C	500 ml	S11M008

\* Radiometer Analytical formulation.

**DEUTSCHER KALIBRIERDIENST DKD**

Calibration laboratory for electrochemical measuring quantities  
 Laboratoire d'étalonnage en mesure électrochimique  
 accredited by the / accrédité par la  
 Akkreditierungsstelle des DKD bei der  
 PHYSIKALISCH-TECHNISCHEN BUNDESANSTALT (PTB)



Calibration certificate  
 Certificat d'étalonnage

Calibration mark  
 Marque d'étalonnage

Object  
 Objet  
**pH calibration standard 10,012**  
**Etalon pH 10,012**

Manufacturer  
 Fabricant  
**HACH LANGE GmbH**  
**D-14163 Berlin**

Nominal value  
 Valeur nominale  
**pH = 10,012 ± 0,010 at 25°C**

Serial / Lot number  
 Numéro de série / lot  
**C01311**

Manufactured for  
 Produit pour  
**Radiometer Analytical SAS**  
**72, rue d'Alsace**  
**69627 Villeurbanne Cedex, France**

Order / Part number  
 Numéro de commande / article  
**S11M007**

Number of pages of the certificate  
 Nombre de pages  
**2**

Date of calibration  
 Date d'étalonnage  
**23.05.2007**

This calibration certificate may not be reproduced other than in full except with the permission of both the Accreditation Body of the DKD and the issuing laboratory. Calibration certificates without signature and seal are not valid.  
 Ce Certificat d'étalonnage ne doit être divulgué que dans sa forme complète et sans modifications. Des extraits ou modifications doivent être autorisés par le Service d'accréditation du DKD et par le laboratoire d'étalonnage ayant établi le certificat. Les certificats d'étalonnage non signés et non estampillés ne sont pas valides.

Seal  
 Cachet  
 Date  
 Date  
 Head of the calibration laboratory  
 Directeur du laboratoire d'étalonnage  
 Person in charge  
 Personne responsable  
 23.05.2007  
 Dr. Helmut Simon  
 Jessica Klinnert



HACH LANGE GmbH  
 D-14163 Berlin  
 ++49 (30) 8 09 86-0  
 ++49 (30) 8 09 86-270

**Certified conductivity standards**

Type	Value	Qty	Part No.
KCl 1D	111.3 mS/cm ±0.5% at 25°C	500 ml	S51M001
KCl 0.1D	12.85 mS/cm ±0.35% at 25°C	500 ml	S51M002
KCl 0.01D	1408 µS/cm ±0.5% at 25°C	500 ml	S51M003
NaCl 0.05%	1015 µS/cm ±0.5% at 25°C	500 ml	S51M004
NaCl 25	25.0 µS/cm ±5% at 25°C	250 ml	S51M013**



\*\* The NaCl 25 certified conductivity standard comes in a 250 ml glass bottle.

**Economical solutions**

For less demanding environments, we offer cost-effective **4-7-10 Series** buffers manufactured to Radiometer Analytical specifications. Radiometer Analytical also produces a range of **molar KCl solutions** for calibrating conductivity cells.

**4-7-10 Series**



Type	Value	Qty	Part No.
pH4.00	pH 4.00 at 25°C	500 ml	S11M012
pH7.00	pH 7.00 at 25°C	500 ml	S11M013
pH10.00	pH 10.00 at 25°C	500 ml	S11M014

**Molar KCl solutions**

Type	Value	Qty	Part No.
KS910	0.1 M KCl (12.88 mS/cm at 25°C)	500 ml	C20C250
KS920	0.01 M KCl (1.413 mS/cm at 25°C)	500 ml	C20C270
KS930	0.001 M KCl (146.9 µS/cm at 25°C)	500 ml	C20C280



## Electrode Maintenance

Proper maintenance of electrodes means accurate and reproducible pH results and longer lifetime. To ensure you get the best possible performance from your electrode, Radiometer Analytical offers a complete kit plus a selection of maintenance and filling solutions.



### GK•ANNEX Electrode Maintenance Kit

The GK•ANNEX Electrode Maintenance Kit (Part no. S91M001) contains all you need for regular maintenance of combined pH, glass and reference electrodes with saturated potassium chloride salt-bridge:

- ✓ Solutions for normal and intensive cleaning,
- ✓ KCl as saturated solution and crystals for refilling,
- ✓ GLP•Logbook to keep note of calibration and maintenance procedures,
- ✓ Utensils and easy-to-follow instructions.

### Maintenance solutions

The maintenance requirements of electrodes vary according to the application. After measurements in aqueous solutions, cleaning with a mild solution such as ○RENOVO•N is sufficient. However, when the solutions contain protein, for example, stronger cleaning agents such as □RENOVO•X or Pepsin in HCl are required. For junctions contaminated with sulphides or AgCl precipitate, Radiometer Analytical's Thiourea Solution is ideal. The operating instructions provide cleaning instructions for each type of electrode.

Type	Description	Part No.
○RENOVO•N	Normal Cleaning Solution, 250 ml	S16M001
□RENOVO•X	Xtra Strong Cleaning Solution, 250 ml	S16M002
KS400	Pepsin in HCl Solution, 250 ml	C20C370
KS410	Thiourea Solution, 250 ml	C20C380

### Filling and storage solutions

Electrode operating instructions contain full information about short and long-term storage of electrodes as well as tips on replenishing the salt-bridge where appropriate. Remember to rinse your electrode thoroughly in distilled water and seal the filling hole before putting it away. The instructions tell you whether to store the sensing element dry or in a particular solution. Follow them carefully as the lifetime of your electrode depends on its correct maintenance and storage.

### Selecting the right filling solution

Reference element	Solutions
Red Rod	KCl•L and KCl•C KS100 and KCl•C
Ag/AgCl (pHC3xxx / REF3xx)	KCl•Ag
Ag/AgCl (XCxxx / XRxxx)	KS120
Calomel	KCl•L and KCl•C KS100 and KCl•C
Hg/Hg <sub>2</sub> SO <sub>4</sub>	KS160

### Filling and storage solutions for combined and reference electrodes

Type	Description	Part No.
KS100	Saturated KCl Solution, 500 ml	C20C300
KCl•L	Saturated Solution of KCl, 100 ml	S21M002
KCl•L-30	Saturated Solution of KCl, 30 ml	S21M010
KCl•C	KCl Crystals, 15 g	S21M001
KCl•Ag	3 M KCl Solution Saturated with AgCl, 100 ml	S21M004
KCl•Ag-30	3 M KCl Solution Saturated with AgCl, 30 ml	S21M011
KS110	KCl Solution 3 M, 500 ml	C20C320
KS120	Saturated KCl + AgCl Solution, 500 ml	C20C310
<del>KS120-30</del>	<del>Saturated KCl + AgCl Solution, 30 ml</del>	<del>S21M009</del>
KS160	Saturated K <sub>2</sub> SO <sub>4</sub> Solution, 500 ml	C20C500
KS160-30	Saturated K <sub>2</sub> SO <sub>4</sub> Solution, 30 ml	S21M007
LiCl-30	1 M LiCl Solution, 30 ml	S21M008

### MeterLab Solution Kit

(Part no. S91M002)

To help you obtain accurate calibrations with your pHC2xx or pHG2xx + REF2xx Red Rod electrodes. It contains: six certified IUPAC pH standards, coloured beakers, beaker dispenser and beaker holder, KCl•L solution and KCl•C crystals.










### pH Meter Solution Kit

(Part no. S91M003)





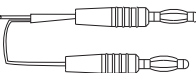


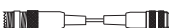

To help you obtain accurate calibrations with your pHC3xxx or REF3xx + pHG3xx electrodes. It contains: three 4-7-10 Series buffers plus KCl•Ag filling solution.

# Accessories

## Connection cables for screw cap electrodes

Instrument socket	Type	Part no.
Banana		CL111 A94L111
<del>P2</del>	<del></del>	<del>CL112 A94L112</del>
DIN		CL113 A94L113
BNC		CL114 A94L114
Type 7		CL116 A94L116
UHF (PL259)		CL119 A94L119
2 x banana		CL120 A94L120
LEMO Ø 5 mm		CL129 A94L129
MAB6		CL136 A94L136

## Plug adapters

Electrode plug	Instrument socket	Part no.
BNC		Type 7 A94P800
Type 7		BNC A94P802
2 x banana		BNC A94P801
BNC		DIN A94P806
Banana		BNC A94P807
<del>Banana Ø 4 mm</del>	<del></del>	<del>Banana US Ø 2 mm X41V001</del>
Banana Ø 4 mm		Type 7 809-147
UHF		MAB6 A94P001
MAB6		UHF+PYE A94P002

## Other accessories

Part no.	Description
X31M012	Conical joints for 12 mm diameter electrodes, polyethylene, set of 4 pcs.
X91M005	Electrode accessory kit consisting of: - 5 protective caps for 7.5 mm diameter electrodes - 5 protective caps for 12 mm diameter electrodes - 5 filling hole clips for 7.5 mm diameter reference electrodes - 5 filling hole clips for 12 mm diameter reference electrodes





## Leading the field in electrochemistry

Radiometer Analytical SAS develops and manufactures an extensive range of electrochemical systems dedicated for routine testing, research and teaching in the laboratory and on the plant.

By supplying instruments, software, sensors and calibration standards, Radiometer Analytical masters the complete measuring chain. Our customers obtain a reliable result at reasonable cost thanks to all-in-one systems that are easy to use and maintain.

The company enjoys a reputation for excellence in the following fields:

**pH, ion and conductivity measurements:** complete systems for reliable measurements in the field and in the lab including a wide choice of instruments, sensors and standards.

**Titration:** workstations customised to individual applications including titrators, sample changers and dedicated software.

**Voltammetry:** all-in systems for electrochemical measurements including potentiostats, impedance meters and powerful software making use of

techniques such as voltammetry, amperometry, coulometry, polarography and EIS.

Radiometer Analytical has been building its expertise for more than seventy years, ever since the company pioneered its very first pH meter in Copenhagen, Denmark. This expertise was strengthened by the acquisition of Tacussel, another leading name in electrochemical instrumentation.

Now based in Lyon, France, Radiometer Analytical today belongs to the Danaher Corporation.

The company is part of the Hach Lange Group whose Calibration Laboratory is accredited for the calibration of pH and conductivity standards by DKD, the German national accreditation authority. Radiometer Analytical SAS is ISO 9001 certified.

Radiometer Analytical products are distributed worldwide through the Hach and Hach Lange sister companies by a network of experienced specialists, offering comprehensive applications and after-sales service.

## Reliable and long-lasting electrodes – the Radiometer Analytical secret

Radiometer Analytical offers a range of more than 300 electrodes - combined pH, glass or reference electrodes, metal electrodes, ion-selective electrodes and conductivity cells - for every application and budget. Electrodes are manufactured on our premises in Villeurbanne, France using a combination of traditional know-how and state-of-the-art technology.

It takes between 2 and 11 days to manufacture a combined pH electrode, depending on the type.

The most spectacular stage of the process is the blowing of the glass bulb from a blob of molten glass heated to 1200°C.

To explain just what goes on behind the scenes when a combined pH electrode is manufactured, Radiometer Analytical has produced an informative illustrated article. Ask for a free copy or download it from our web site: [www.hach-lange.com](http://www.hach-lange.com).



Preparing the stem for dipping



Dipping in molten glass



Blowing the glass bulb

Data subject to change without notice.

Distributed by



Australian Head Office

PO BOX 465  
Blackburn, VIC 3130

Free Call: 1800 723 405

Ph: (03) 9845 0300

Fax: + 61 3 9845 0350

e-mail: [sales@pacificlab.com.au](mailto:sales@pacificlab.com.au)

Connect with us:



/PacificLab



@PacificLabP

