

# Axygen® Automation-compatible PCR Microplate

Unique design prevents plate warping and minimizes evaporation during thermal cycling

- Better sample traceability
- Equivalent performance
- Alternative product


Warping of PCR microplates through thermal cycling increases sample evaporation and makes them no longer automation-compatible. As a result, Corning invented a unique, one material automation-compatible PCR microplate that sustains flatness post-thermal cycling, has raised-angled rims for improved sealing, and has grid lines for better sample traceability. This microplate is compatible with major automated thermal cyclers, liquid handlers, and sequencers.

## Features

- ▶ A rigid frame with a small split **to remain flat** during thermal cycling without distorting its shape.
- ▶ Specially designed raised rims around the wells provide tight sealing to **minimize evaporation**
- ▶ Grid lines, clear wells, and alphanumeric lettering for improved **sample traceability**
- ▶ Compatible with most major thermal cyclers, sequencers, and automated liquid handlers
- ▶ Ultra-thin walls provide fast and consistent heat transfer during PCR.
- ▶ Nonpyrogenic
- ▶ RNase-/DNase-free
- ▶ PCR inhibition and human DNA-free

**AXYGEN**

A Corning Brand

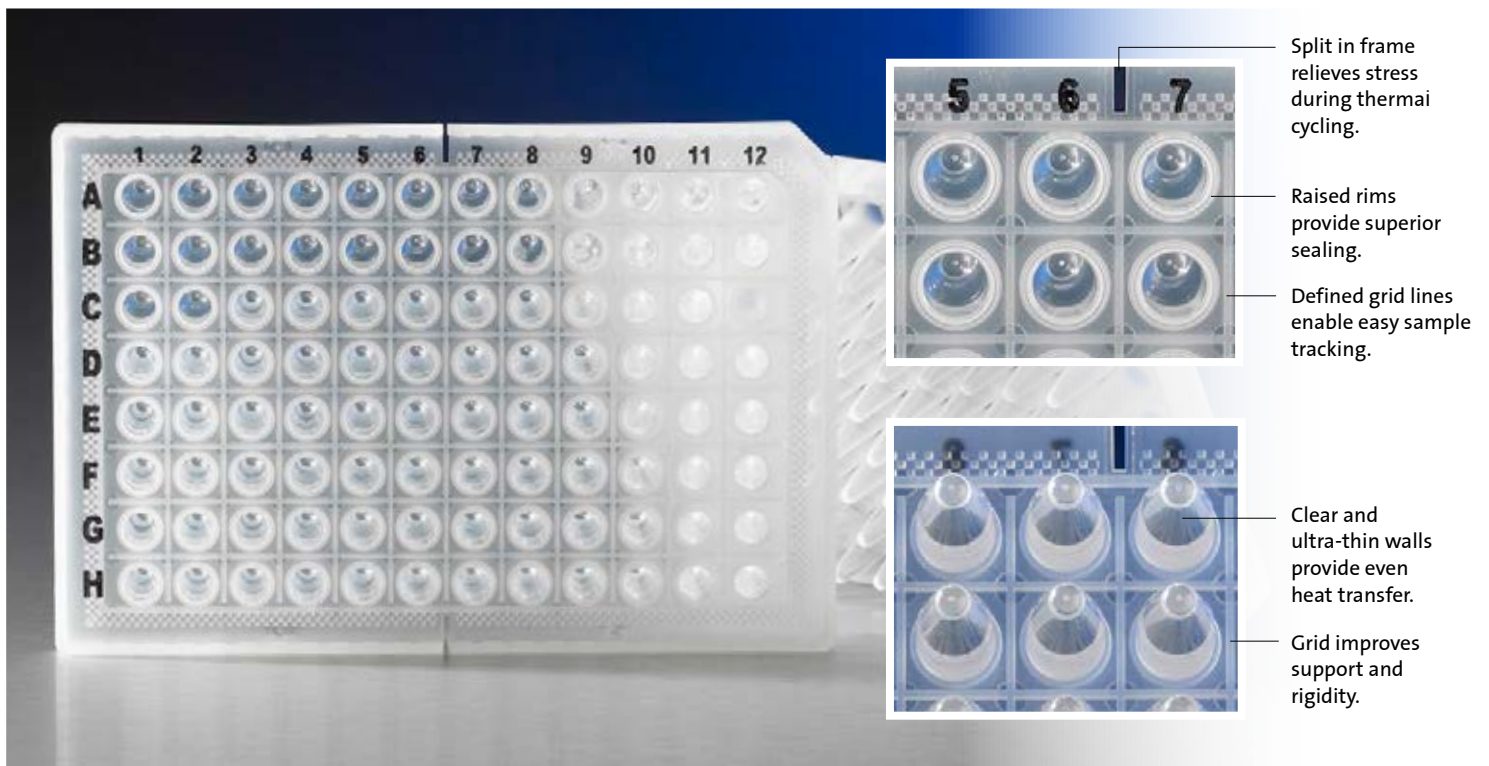

## Ordering Information

### Axygen® Automation-compatible PCR Microplate

Cat. No.	Description	Volume (µL)	Qty/Pk	Qty/Cs
PCR-96-HS-AC-C	Axygen 96-well half skirt automation-compatible, PCR microplate, clear, nonsterile	200	10	50

### Accessories

Cat. No.	Description	Qty/Pk	Qty/Cs
AM-96-PCR-RD	Sealing mat for 96-well PCR plate	10	50
UC-500	UltraClear pressure sensitive film	100	500
PCR-02-FCP-C	0.2 mL clear, thin wall PCR 8-strip flat caps with writing surface, nonsterile	125	1,250
PCR-02CP-C	0.2 mL clear, thin wall PCR 8-strip domed caps, nonsterile	125	1,250
PCR-2CP-RT-C	0.2 mL ultra-clear, thin wall PCR 8-strip caps for real time PCR, nonsterile	125	1,250
PCR-12-FCP-C	0.2 mL clear, thin wall PCR 12-strip flat caps, nonsterile	80	800

### Compatible with most major real-time PCR thermal cyclers, including:

- ▶ Applied Biosystems 7000, 7300, 7500, 7700, 7900HT, ProFlex™, StepOne™, StepOnePlus™, ViiA™ 7, QuantStudio™ 12K Flex
- ▶ BioRad/MJ Research iCycler®, MyiQ™, iQ5™, Chromo4™, DNA Engine Opticon® 2, CFX96 Touch™, CFX Connect™, MiniOpticon™
- ▶ Eppendorf Mastercycler® ep realplex

### Compatible with most major standard PCR thermal cyclers, including:

- ▶ Applied Biosystems Veriti®, 9600, 9700, 2720
- ▶ BioRad/MJ Research MyCycler™, T100™, C1000 Touch™, S1000™, DNA Engine Tetrad® 2, PTC-100®, PTC-200®, PTC-225 Tetrad, DNA Engine Dyad®, Dyad Disciple™
- ▶ Eppendorf Mastercycler series: ep, pro, nexus, nexus eco, nexus gradient, nexus gradient eco, nexus GSX1, nexus GSX1e, nexus SX1, nexus SX1e

### Compatible with the following Applied Biosystems sequencers:

- ▶ 3100, 3130, 3700, 3730/3730x

### Compatible with most major automated liquid handlers, including:

- ▶ Beckman Coulter, Hamilton, Tecan, Agilent, Perkin-Elmer

For more specific information on claims, visit the Certificates page at [www.corning.com/lifesciences](http://www.corning.com/lifesciences).


**Warranty/Disclaimer:** Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

**CORNING** | **FALCON** | **AXYGEN** | **GOSELIN** | **PYREX**



Australian Head Office  
PO BOX 465  
Blackburn, VIC 3130  
Free Call: 1800 723 405  
Ph: (03) 9845 0300  
Fax: + 61 3 9845 0350  
e-mail: [sales@pacifilab.com.au](mailto:sales@pacifilab.com.au)

Connect with us:

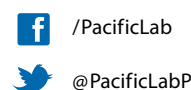

[www.pacifilab.com.au](http://www.pacifilab.com.au)