


EndoSwab | Endoscope Sample Collection



Hygiena EndoSwab enables cleaning verification of scopes or other reusable instruments with hard to reach internal channels. The device features an absorbent foam tip attached to a long flexible wand. After sampling, the swab tip is removed and used with Hygiena ATP tests to detect any residual contamination left behind after cleaning, ensuring scopes and other reusable instruments are ready for disinfection and cleaning procedures were effective.

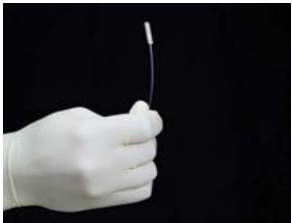


Benefits

- Results in 15 seconds
- Various diameter absorbent sponge tip for maximum sample recovery
- Length ensures recovery of a representative sample for most endoscope lengths
- Clinically clean packaging
- Individually bagged for single use

Compatibility

- Use with Hygiena UltraSnap ATP tests or SuperSnap high sensitivity ATP tests, depending on desired sensitivity.

EndoSwab Procedure

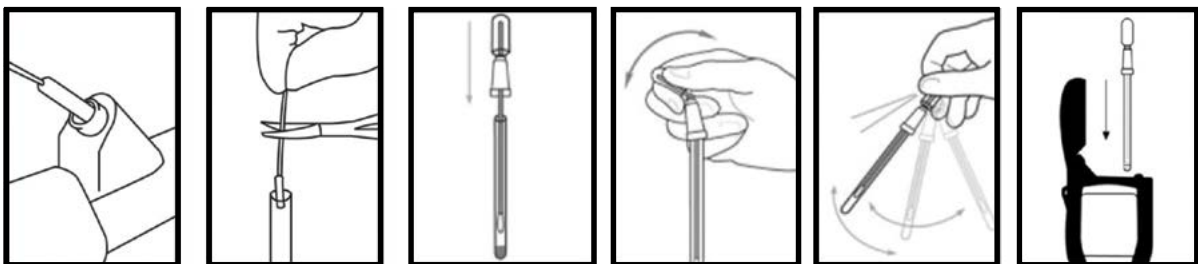

Part Number	Description		Quantity
HG-ENDO1MM	Endoswab – Swab for Endoscope Cleaning Verification, 1mm diameter	Shop	100
HG-ENDO2MM	Endoswab – Swab for Endoscope Cleaning Verification, 2mm diameter	Shop	100
HG-ENDOSWAB	Endoswab – Swab for Endoscope Cleaning Verification, 3mm diameter	Shop	100
HG-ENDO4MM	Endoswab – Swab for Endoscope Cleaning Verification, 4mm diameter	Shop	100